




Sequential Integrator Class 5



SAP / Monitoring

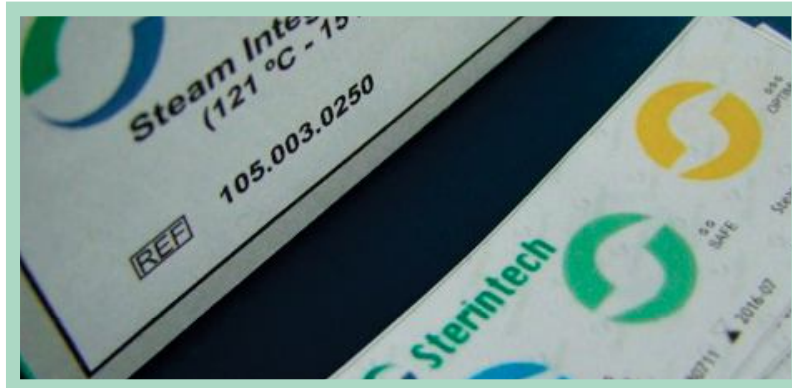

Integrator with a sequence

REF: 105.003.0250 Sequential Integrator 134°C-3,5 min Unit: 250 pieces
REF: 105.018.0250 Sequential Integrator 134°C- 18 min Unit: 250 pieces

Integrators can be used in addition to the batch control Helix which will allow you to release load. Especially when OR staff wants to be re-affirmed about the sterilization conditions inside trays or containers.

The Steritech™ Integrators are offering you more information about the steam sterilization process than those with only one integrator printed. Due to that the three integrators are calibrated to different time/temperature combination it is providing information about for instance over-heating.

Characteristics:

- Shelf life: 5 years after production
- Packed in an easy dispensing carton
- Steam Integrator conform to ISO 11140-1 : Class 5
- Clear colour transition into black of the integrator
- Chemical Integrators covered with lacquer
- Non-Toxic / Lead Free

