




# Bowie Dick Test Pack



## SAP / Monitoring

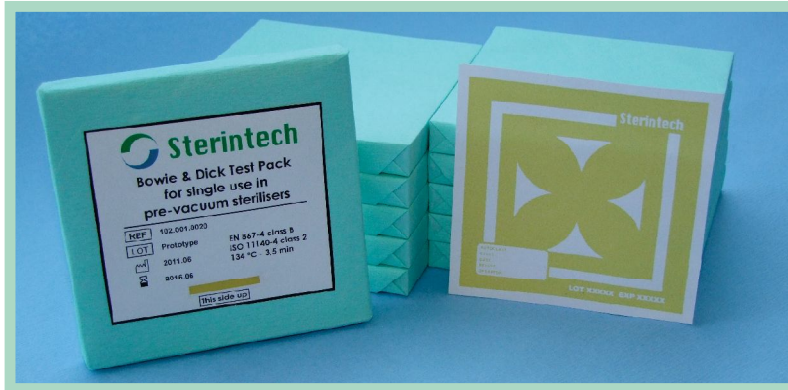

## Once a day ..... since 1963

REF: 102.001.0020 Unit: 20 pieces

Used for daily monitoring of pre-vacuum steam sterilizers in accordance to EN 285 requirements. Based upon the Bowie Dick test result the steam sterilizer is released for usage during the day.

### Characteristics:

- Fulfil requirements of EN 285 and ISO 17665
- Conform to ISO 11140-1 and ISO 11140-4 : Class 2'
- Independently tested to show performance similar to original pack
- Clear yellow – black colour change
- Easy to interpret
- Easy to document
- Non-Toxic / Lead Free
- Dimensions: 120 x 120 x 15 mm

The Bowie Dick Test Pack Detect failures like:

- Air leaks,
- Insufficient Vacuum,
- Inadequate Steam Penetration.

